

Dental Handpiece Maintenance Unit

LUBRICARE

Operation Instructions

Thank you for purchasing LUBRICARE, a handpiece maintenance unit of Bien-Air Dental SA.
LUBRICARE lubricates and internally cleans dental handpieces easily and effectively.
For optimum safety and performance, read this manual thoroughly before using the unit and pay close attention to the warnings and notes.
Keep this manual in a handy place for ready reference.

Table of Contents

	Page
PREVENT ACCIDENTS.....	2
ATTENTION.....	3
1. Warning and Prohibitions.....	4
2. Technical Specifications	5
3. Parts Identification	8
4. Set Up.....	11
5. Operation	16
6. Maintenance and Replacement Parts	25
7. Maintenance and Inspection	30
8. Troubleshooting.....	31
9. Service Contacts	32
10. Appendix – Electromagnetic Declaration.....	33

ATTENTION CUSTOMERS

Do not fail to receive clear instructions concerning the various ways to use this equipment as described in this accompanying Operator's Manual.
Fill out and sign the warranty and give a copy to the dealer from whom you purchased the equipment.

ATTENTION DEALERS

Do not fail to give clear instructions concerning the various ways to use this equipment as described in this accompanying Operator's Manual.
After instructing the customer in the operation of the equipment, have the customer fill out and sign the warranty. Then fill in your own section of the warranty and give a copy to the customer. Do not fail to send the manufacturer's copy to Bien-Air Dental SA regional office.

PREVENT ACCIDENTS

Most operation and maintenance problems result from insufficient attention being paid to basic safety precautions and not being able to foresee the possibilities of accidents. Problems and accidents are best avoided by foreseeing the possibility of danger and operating the unit in accordance with the manufacturer's recommendations. First thoroughly read all precautions and instructions pertaining to safety and accident prevention; then, operate the equipment with the utmost caution to prevent either damaging the equipment itself or causing bodily injury.

Note the meaning of the following symbols and expressions:



This warns the user of danger of death, serious bodily injury or total equipment damage and failure or fire.



This identifies methods which must not be used or purposes which the instrument is not suited for.



This alerts the user to the risk of light to medium injury or equipment damage.



This alerts the user of important points concerning operation.

The user (e.g., healthcare facility, clinic, hospital etc.) is responsible for the management, maintenance and use of medical devices. Also this equipment must not be used by anyone except legally qualified dentist, doctor, or other legally qualified professional.

Do not use this equipment for anything other than its specified dental treatment purpose.

Rx Only

Prescription Device

CAUTION: Federal law restricts this device to sale by or on the order of a dentist. (U.S.A.)

ATTENTION

- Bien-Air Dental SA will not be responsible for accidents, instrument damage, or bodily injury resulting from:
 1. Repairs made by personnel not authorized by Bien-Air Dental SA.
 2. Any changes, modifications, or alterations of its products.
 3. The use of products or instruments made by other manufacturers, except for those procured by Bien-Air Dental SA.
 4. Maintenance or repairs using parts or components other than those specified by Bien-Air Dental SA. and other than in their original condition.
 5. Operating the instrument in ways other than the operating procedures described in this manual or resulting from the safety precautions and warnings in this manual not being observed.
 6. Workplace conditions and environment or installation conditions which do not conform to those stated in this manual such as improper electrical power supply.
 7. Fires, earthquakes, floods, lightning, natural disasters, or acts of God.
- The useful life of LUBRICARE is 5 years (based on self-certification) from the date of shipment provided it is regularly and properly inspected and maintained.
- Bien-Air Dental SA will supply replacement parts and be able to repair the product for a period of 10 years after the manufacture of the product has been discontinued.

1. Warnings and Prohibitions

WARNING

- *No modification of this equipment is allowed. The customer is not allowed to repair this equipment himself.*

PROHIBITION

- *LUBRICARE cannot be used for air bearing handpieces (Astron series).*
- *Electromagnetic wave interferences caused by cellular phones, transceivers, remote controls and similar transmission devices could cause the equipment to operate randomly. All devices which transmit electromagnetic waves located near the work area should be turned off.*

2. Technical Specifications

Specifications

Name	LUBRICARE
Model	HIM-1
Type	BA
Rating	AC 100 V to 240 V 50/60 Hz
Power Consumption	25 VA
Fuse	250 V 2 A Slow Blow and High breaking Type $\varnothing 5 \times 20$ mm
Input Air Pressure	1.0 MPa max.
Air Pressure	0.3 MPa to 0.5 MPa
Class	Class I
Isolation from the Supply Mains	Unplug the power cord from supply mains
Recommended Air Pressure	0.35 MPa
Air Flow Rate	40 NL/min to 60 NL/min
Weight	Approx. 10 kg
Size (including regulator)	Width 300 mm \times Height 370 mm \times Depth 300 mm

Product Description

Used to maintain optimum performance and prolong working life of dental handpieces.
Delivers oil and air automatically to handpiece.
Used after dental treatment and before autoclaving.

Operating Principle

Physical methods used to accomplish its intended use:

Air Pressure: 0.3 MPa to 0.5 MPa
Air Flow Rate: 40 NL/min to 60 NL/min
Maintenance oil spray can

Mechanisms by which it works:

Maintenance spray can is operated by pressurized air. Sprayed air delivers oil and lubricates the inside of the handpieces. Also, excess oil is removed from the handpieces.

Intended User

- a) Engaged person (peoples) in dental clinic
- b) Language Understanding: English or languages offered in the instruction for use. Understanding of attention and warning marks.
- c) Experience: Not relevant

Operating Environments

Dental clinic, hospital
Temperature: +10°C to +40°C (+50°F to +104°F)
Humidity: 30% to 75% (without condensation)
Atmospheric Pressure: 70 kPa to 106 kPa

Transport and Storage Environments

Temperature: -10°C to +50°C (+14°F to +122°F)
Humidity: 30% to 75% (without condensation)
Atmospheric Pressure: 50 kPa to 106 kPa

Disposal

The package should be recycled. Metal parts of the equipment are disposed as scrap metal. Synthetic materials, electrical components, and printed circuit boards are disposed as electrical scrap. Material must be disposed according to the relevant national legal regulations. Consult specialized disposal companies for this purpose. Please inquire of the local city / community administrations concerning local disposal companies.

Meaning of Symbols

* Some symbols may not be used.



Manufacturer



Date of manufacture



Serial number



Attention, consult accompany documents



Refer to instruction manual



Alternating current

INPUT

Rated input amperes

U

Rated supply voltage

f

Supply frequency



Non-ionizing electromagnetic radiation



Main switch
Off On



Stand-by
If the main switch is on, green power lamp next this mark lights up.



General warning sign



Warning; Flammable material



No open flame; Fire, open ignition source and smoking prohibited



GS1 DataMatrix



This way up



Fragile; handle with care



Keep away from rain



Temperature limitation



Humidity limitation



Atmospheric pressure limitation

Rx Only

Prescription Device
CAUTION: Federal law restricts this device to sale by or on the order of a dentist. (U.S.A.)



cTUVus certification mark (U.S.A. and Canada)



CE marking
Conforms with the European Directive, 93/42/EEC and 2011/65/EU.



Authorized representative in the European Community



Marking of electrical equipment in accordance with the European directive 2012/19/EU (WEEE)





C-Tick mark
Conforms with applicable electromagnetic compatibility (EMC) requirements (Australia and New Zealand)

N13700

C-Tick supplier code number


3. Parts Identification

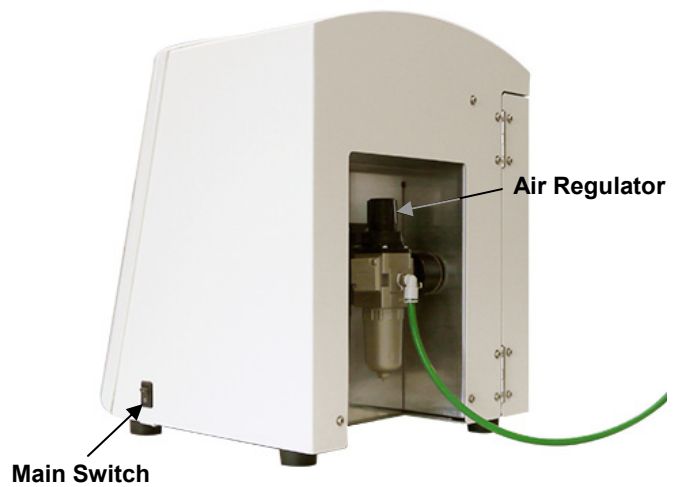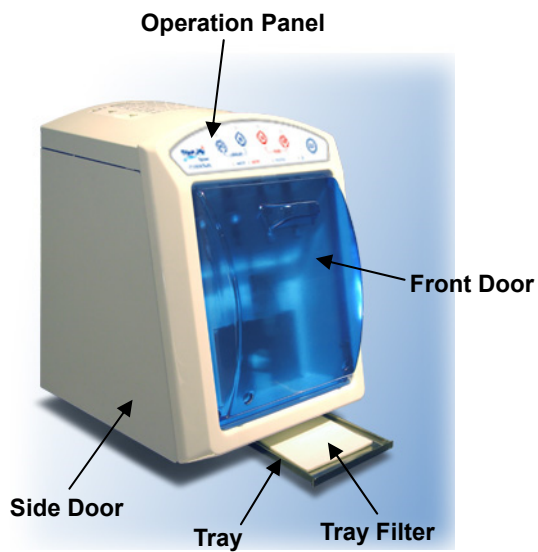

Power Cord Connection

Fuse Holder


Rating Label

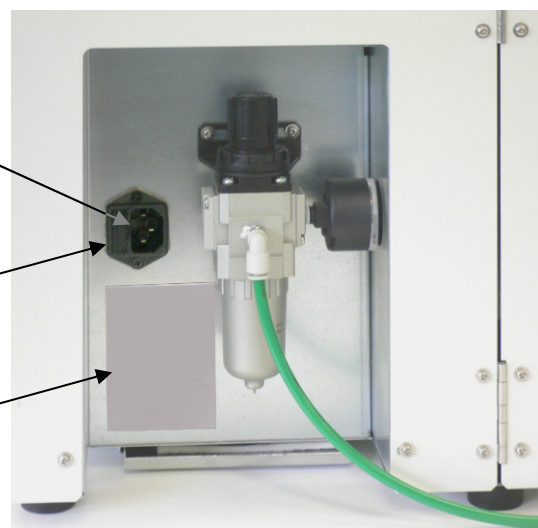

Selector Switches

Selector Lamps (green)



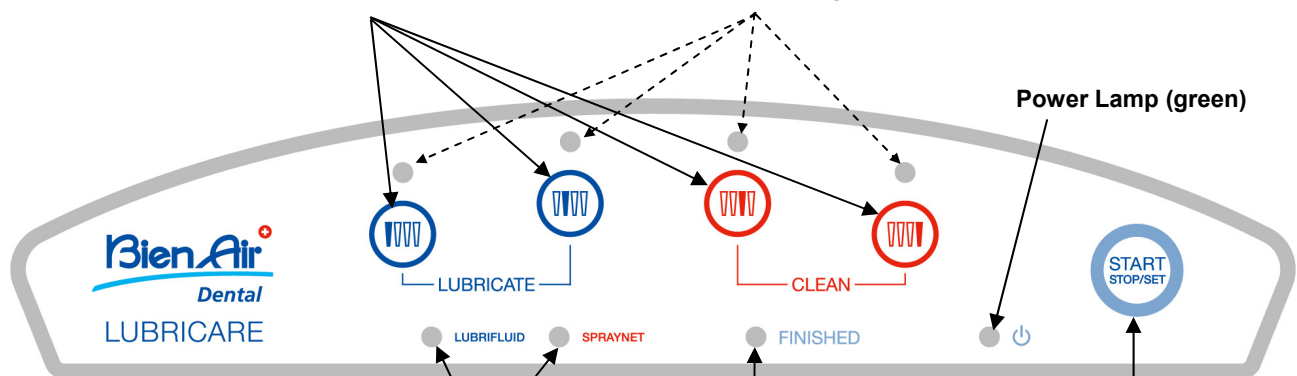


Power Lamp (green)




Spray Can Empty Lamps (orange)

Finished Lamp (green)

Start, Stop and Set Switch

	<div> <div>WARNING</div> <div>DANGER</div> </div>
	<div> <div> <ul style="list-style-type: none"> – Spray cans will explode if they get too hot. – Do not use this device near an open flame or in a spot where it might heat up. – Use only the sprays approved by Bien-Air Dental SA. </div> <div> PREVENT EXPLOSION AND IGNITION OF SPRAY CANS! </div> <div>  </div> <div> <p>Les bombes aérosols peuvent exploser si elles sont exposées à une source de chaleur. Ne pas utiliser ce dispositif à proximité directe d'une flamme ou d'une source de chaleur. Utiliser uniquement les bombes aérosol de Bien-Air Dental SA. EVITER EXPLOSION ET ECHAUFFEMENT DES BOMBES AEROSOLS</p> </div> </div>
	<div> <div> <ul style="list-style-type: none"> – Ventilate the room while using the Lubricare. </div> <div> BREATHING CONTAMINATED AIR IS HARMFUL. </div> <div>  </div> <div> <p>Aérer les locaux durant l'utilisation du Lubricare. RESPIRER L'AIR CONTAMINE EST NOCIF.</p> </div> <hr/> <div> <ul style="list-style-type: none"> – Read the user's manual and follow the instructions carefully. </div> <div> <p>Lire le mode d'emploi et suivre rigoureusement les instructions.</p> </div> </div>

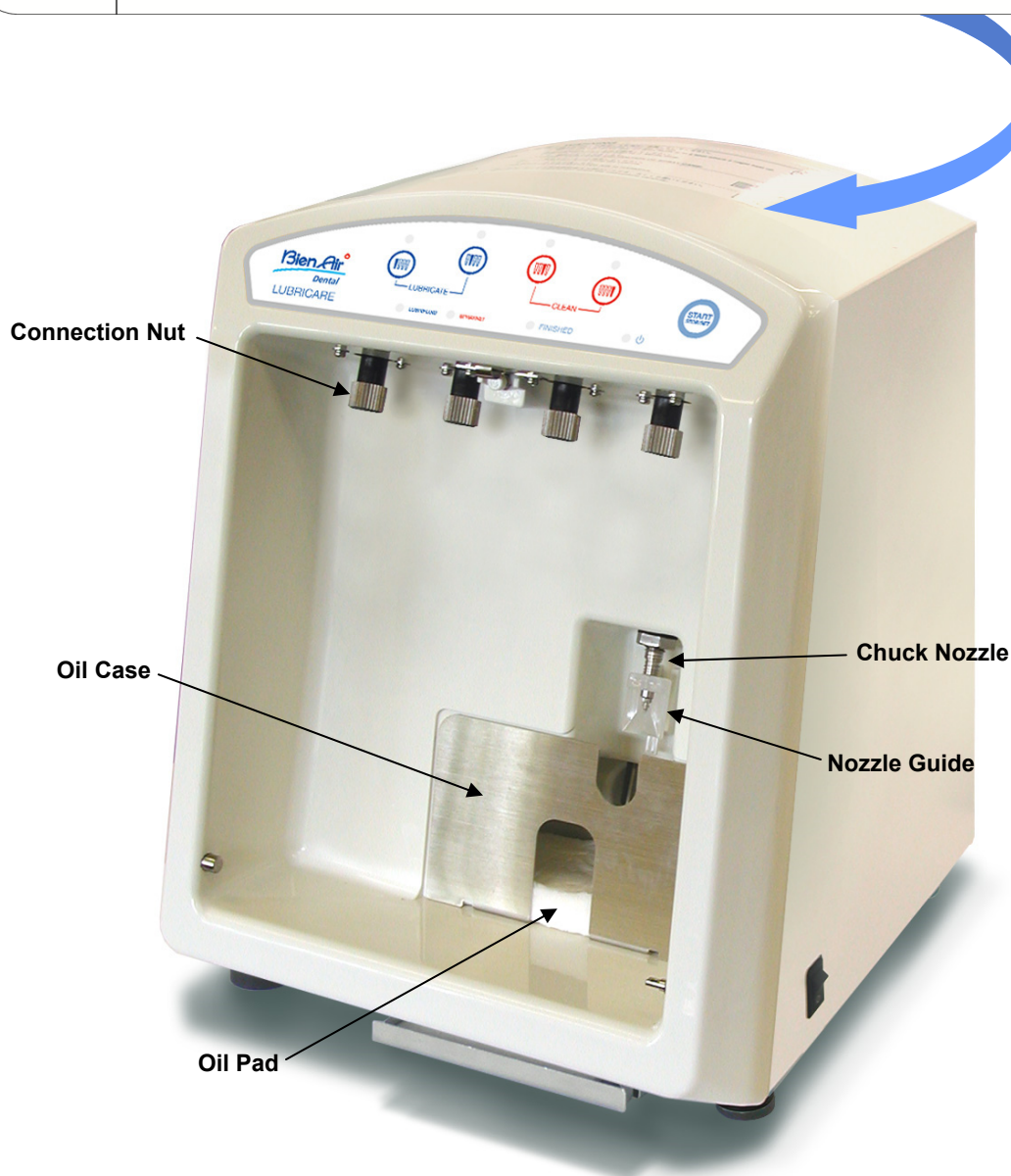




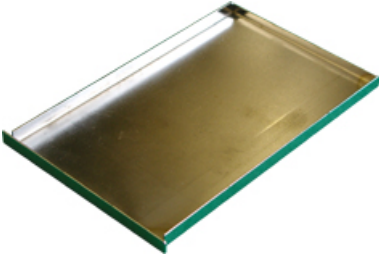

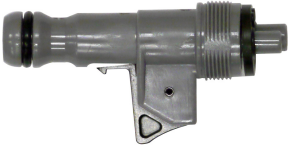



Components

Main Unit



Accessories

Oil Absorption Pads (5)	Oil Absorption Case (1)	Power Cord (1)
		<div>(CE · 230 V) </div> <div>(US · 120 V) </div>
Tray Filter (1)	Tray (1)	Spray Can Stand (2)
		
TR Coupling (2)	Air Tube (1)	
		

4. Set Up

- * Direct inquiries concerning types of couplings, sprays, and spray can stands to your local dealer or Bien-Air Dental SA.

WARNING

- **Explosion Hazard:** Do not use near open flames or other ignition sources.
- Use only spray provided by Bien-Air Dental SA.
- **Health Hazard:** The area must be well ventilated. Breathing fumes could damage your health.
- The equipment must be properly grounded.


CAUTION

Set the unit on a level and stable surface.

Usage Note

- ◆ Leave at least 5 cm of free, open space around the unit.

(1) Air Connection

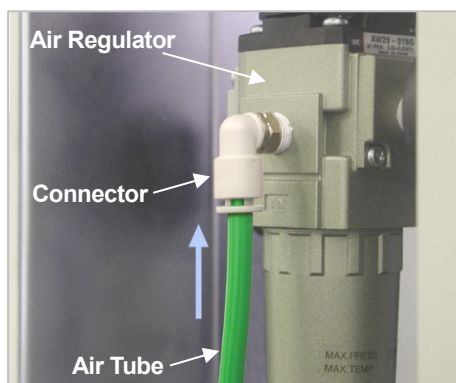

1. Connect the green air tube provided to the air regulator and the main air supply.
Open the valve for the main air supply.

Usage Note

- ◆ The main air supply must be clean and oil-free.
- ◆ Make sure the tube is securely connected.



2. Set the air regulator dial between 0.3 and 0.5 MPa.
Pull the dial up to unlock it. Push it back down after setting the pressure.

- * Recommended Air Pressure is 0.35 MPa.


Usage Note

- ◆ Less than 0.3 MPa will result in poor performance.
- ◆ More than 0.5 MPa could result in damage.

CAUTION

Depending on the type of handpiece and coupling, there could be a release of oil mist. In this case, lower the pressure to 0.3 MPa.

(2) Power Cord

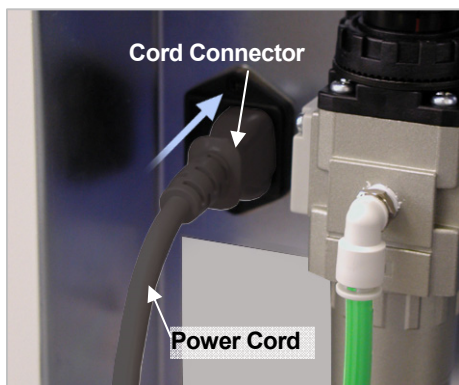


Connect one end of the power cord to the main unit and plug the other end into a standard receptacle.

WARNING

- *Use the right cord!! Since power supply ratings differ depending on the country, two power supply cords are provided. Use the one suited for your local power supply; using the wrong cord could cause a fire.*
- *Make sure the ground is properly connected.*
- *Do not place objects around the power input plug at any time so that the equipment can be cut off from the supply mains when necessary.*
- *Do not place this equipment where the main switch cannot be operated.*

(3) Coupling

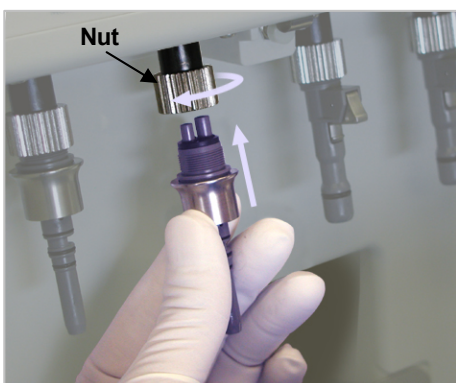


Fit the coupling onto the nozzle and tighten up its nut.


Usage Note

- ♦ *Make sure the nut is properly tightened. Otherwise oil and air will escape and the maintenance will not be properly performed.*

(4) Oil Pad

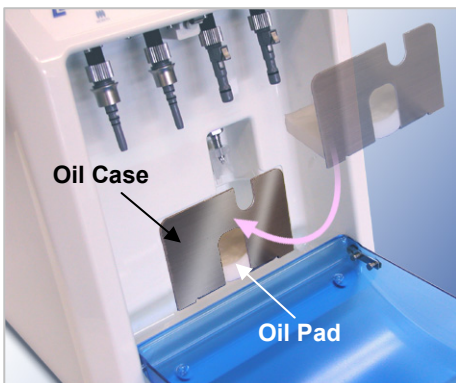

Put an oil pad in its case and install the case.

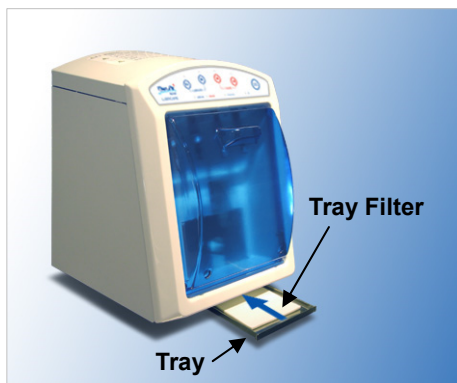

(5) Tray

Put a filter in the tray and slide it into place.


Usage Note

- ♦ *Slide the tray all the way in; otherwise, a film of oil will be deposited on the unit.*


(6) Spray Cans




1. Open the side door.
Press the spot indicated in the photo to open the door.

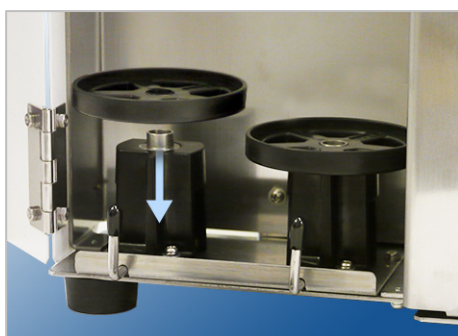

2. Set the stands for the cans in place as shown in the photo.

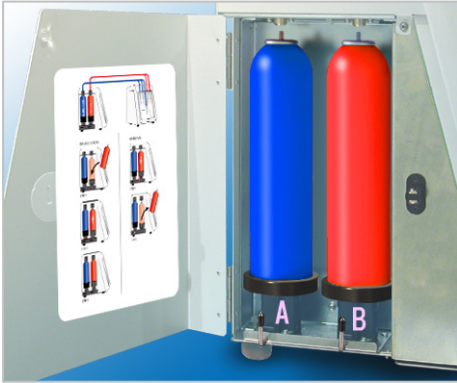

3. Can A is for the chuck maintenance and handpiece nozzles 1 and 2. Can B is for nozzles 3 and 4.

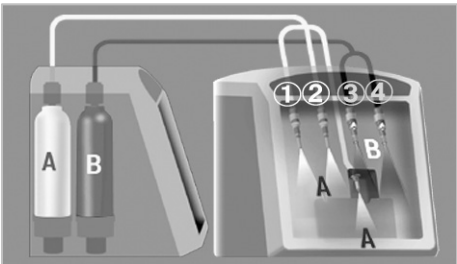

4. Push the lock levers down to a horizontal position to secure the cans.

⚠ CAUTION

Oil will leak or spray out if the cans are not properly installed.

Usage Note


- ♦ *Make sure the lock levers are all the way down. Otherwise, the oil will not come out. (The empty lamp will start blinking.)*
- ♦ *Pull and push the cans to make sure they are secure and do not wobble.*




5. Close the side door.

(7) Spray Can Removal

1. Pull the lock lever straight up to lower the stand.

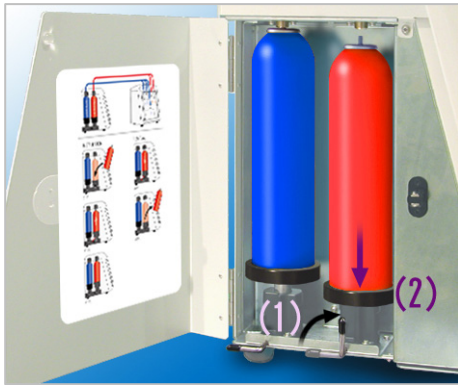

2. Take the can out.
If the can's nozzle is stuck in place, simply pull it down slightly and take it out.

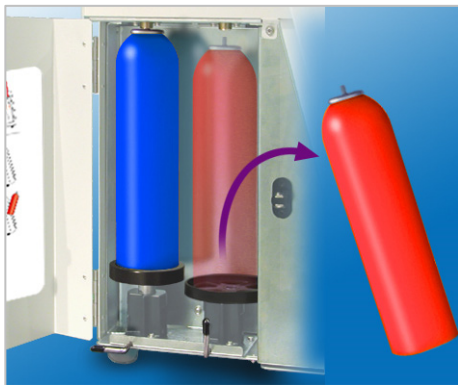


5. Operation

Handpiece Autoclave Sterilization

WARNING

To prevent the spread of grave, life-threatening diseases like HIV and hepatitis B, autoclave handpieces after performing regular maintenance using LUBRICARE.

<Take Off>

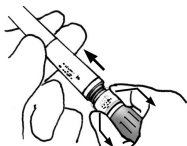

- (1) Take the handpiece off its tube.
(Remove the burr as well.)

<Clean>




- (2) Clean the head with a brush and running water.



- (3) Wipe the handpiece with gauze dampened with Ethanol for Disinfection (Ethanol 70 vol% to 80 vol%).

LUBRICARE Procedure

< Clean and Lubricate >

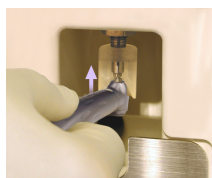

Perform chuck maintenance for Morita ball bearing handpieces.



Perform the maintenance procedure.

<Clean>



Remove excess oil by wiping with Ethanol for Disinfection (Ethanol 70 vol% to 80 vol%).

<Autoclave>

Autoclave in a sterilization pouch.

Usage Note

- ◆ Some oil may remain inside the handpiece even after autoclaving. Refer to the user instructions for each handpiece for proper storage. (Storing the handpiece in an upright position is recommended.)

- * If LUBRICARE hasn't been used for a while, make sure it operates properly and safely before using it.

WARNING

- *Avoid the risks of electrical shock, equipment damage and fire during an electrical storm: Turn LUBRICARE off and do not touch it or its cord.*
- *Wear surgical gloves to operate and clean LUBRICARE.*

<Set Up>

- Carefully clean off the outside of the handpiece before performing maintenance.
- Remove the burr.


Chuck Cleaning and Lubrication

Turn the main switch on.
(Green power lamp lights up.)



Open the front door.




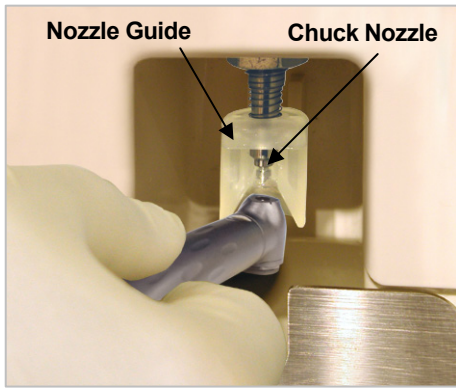


Line up the front end of the head with the nozzle guide and set the chuck in place.

⚠ CAUTION

Make sure the head is flat against the nozzle guide. If it is tilted, the nozzle might be bent or oil might spray out into your eyes.

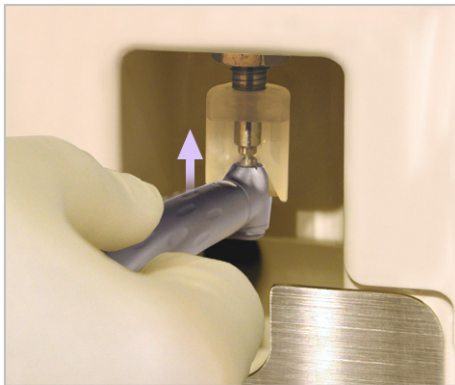

Lift the handpiece up slightly to begin cleaning and lubricating the chuck.

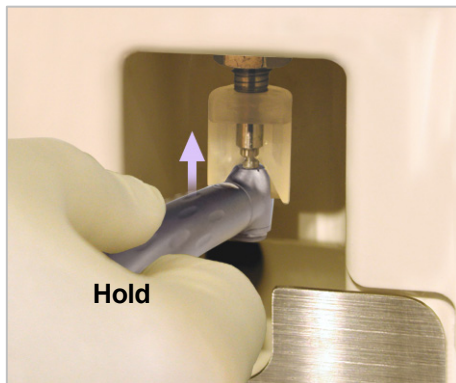

Hold it in place for about 5 seconds.

A beep sounds when the procedure is finished.


* A series of beeps signals an error. In this case, try again.

⚠ CAUTION

Oil will spray out of the nozzle if the handpiece is removed too soon (before the beep).

Usage Note

- ◆ *Chuck cleaning and lubrication is only for Bien-Air handpieces.*
- ◆ *If the nozzle guide comes off, line up the tab on it with the groove and slide it back on.*

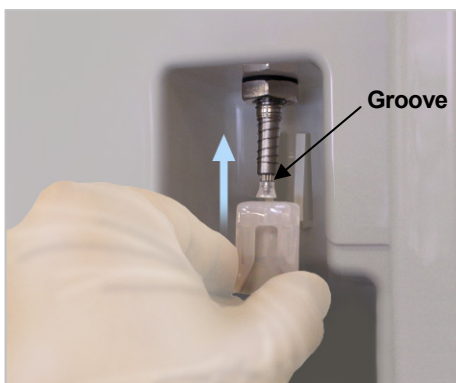

Handpiece Maintenance



- * It takes about 25 seconds to clean and lubricate one handpiece.

Turn on the main switch.
(Green power lamp lights up.)




Open the front door.

- * For micromotor attachments (straight, contra etc.), some should have a burr inserted and some should not. Refer to the user instructions for the micromotor attachments.

CAUTION

- *For some handpieces, an oil mist may be released if a burr is not inserted.*
- *Depending on the type of handpiece and coupling, there could be a release of oil mist. In this case, lower the pressure to 0.3 MPa.*

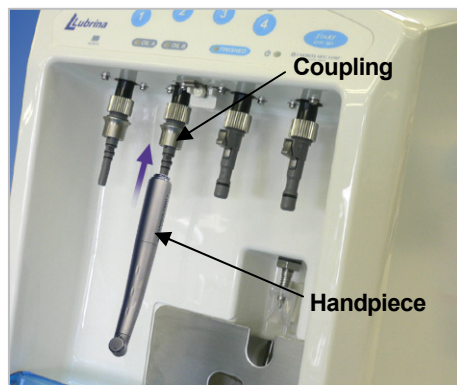

Attach a handpiece to its coupling with the chuck facing back.
(See photo.)

Spray can A is for nozzles 1 and 2. Spray can B is for 3 and 4.

CAUTION


*Do not damage the O-ring on the coupling when attaching a handpiece.
This could result in oil spraying out or poor performance.*

WARNING

- *Do not select a coupling that does not have a handpiece attached to it. Otherwise, you may inhale sprayed oil mist, or it may splash into your eyes.*
- *Select the most appropriate lubrication mode for each handpiece. Otherwise, too much oil will remain inside the handpiece and it will spray out when using the handpiece.*
- *Use this equipment with proper air pressure (0.3 MPa to 0.5 MPa). Otherwise, too much oil will remain inside the handpiece and it will spray out when using the handpiece.*

Usage Note

- ♦ *Put the handpiece all the way on until it clicks into place. A handpiece could come off during the maintenance if it is not properly connected.*
- ♦ *If the chuck is facing the front door, the inside of the unit will tend to get dirtier.*

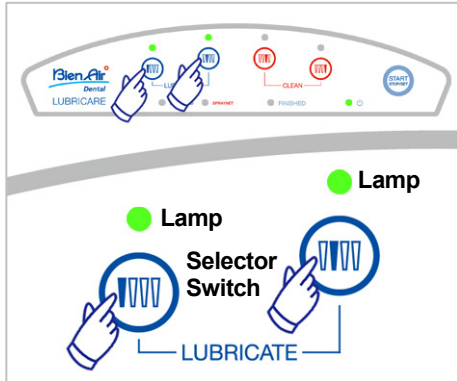


Press the selector switch corresponding to the handpieces which are connected.
Selector lamps light up.
(If you press the wrong one, just press it again to turn it off.)



Close the front door.

Usage Note

- ♦ *LUBRICARE will not start if the front door is not closed.*

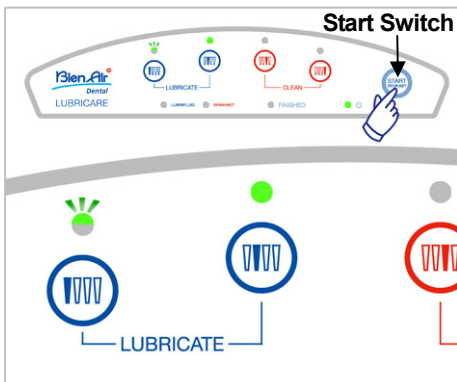

Press the Start Switch.
The handpiece lamps blink while maintenance is performed.

- * Press the Start Switch again if the procedure was stopped before finishing.

CAUTION


Do not select a coupling that does not have a handpiece attached to it. This would result in oil spraying out when the Start Switch is pressed.

Usage Note

- ♦ *Never remove the spray can during maintenance.*
- ♦ *Opening the front door will immediately stop the maintenance procedure.*

- * If a can runs out during maintenance, the empty lamp will start flashing, and the procedure will stop. Replace the can and press the Start Switch to repeat the procedure from the beginning.

When the maintenance for all the handpieces is completed, their lamps will go out, a beep will sound and the Finished Lamp will light up.

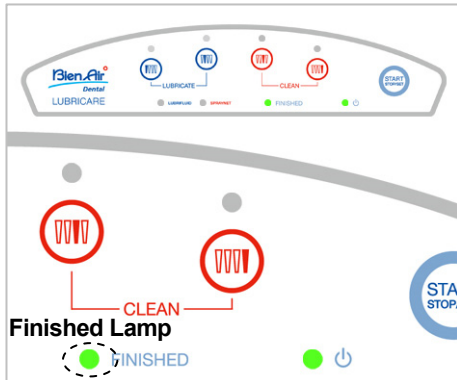

Open the front door and remove the handpieces.
Hold the handpiece with one hand remove it from coupling.

Usage Note

- ♦ *Use Ethanol for Disinfection (Ethanol 70 vol% to 80 vol%) to wipe oil off the handpiece.*




For micromotor attachments, hold the attachment with one hand and press the release button with the other.



Take out the handpieces and close the door.

⚠ CAUTION

Take care not to cut your fingers on the front door or other parts.

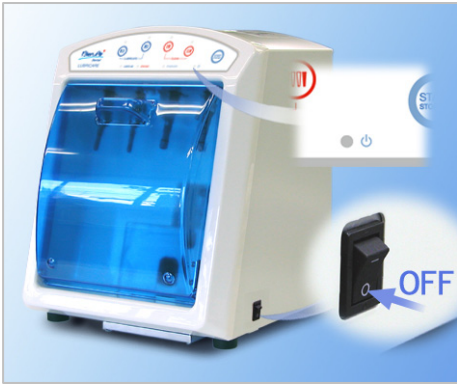

Turn the main switch off after use.
The power lamp will go out.

⚠ CAUTION

- *Do not fail to turn off the main switch after use.*
- *If the unit will not be used for a while, unplug it and close the main air valve.*

Change Lubrication Mode (time)

* LUBRICARE can be used with sprays other than those made by Bien-Air Dental SA. To do this the stand for the spray can and the coupling must be replaced.
Change the lubrication time for other sprays in the following way.

There are three lubrication modes: the short mode (with less oil), the long mode (with more oil), and the long blow-out mode (to remove excess oil more thoroughly).

Short Mode

* For less oil

The selector switch lamp blinks at 0.3 second intervals.
Lubrication time for one handpiece is about 25 seconds.

Long Mode

* For more oil

The selector switch lamp blinks at 1.0 second intervals.
Lubrication time for one handpiece is about 25 seconds.


Long Blow-out Mode

* To remove excess oil more thoroughly


The selector switch lamp blinks at 2.0 second intervals.
Lubrication time for one handpiece is about 120 seconds.
This mode blows out excess oil for a longer time after the short mode completes the oiling procedure.

Usage Note


♦ *For some types of spray cans, the long setting will supply too much oil and it will spray out of the head.*

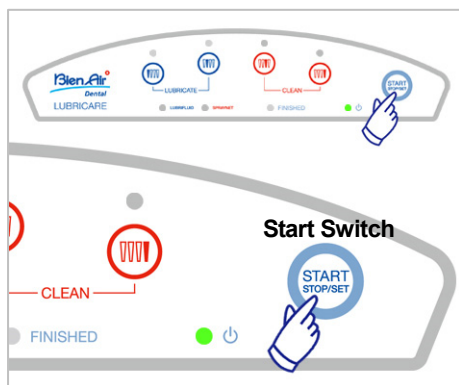

1. Turn the main switch on.

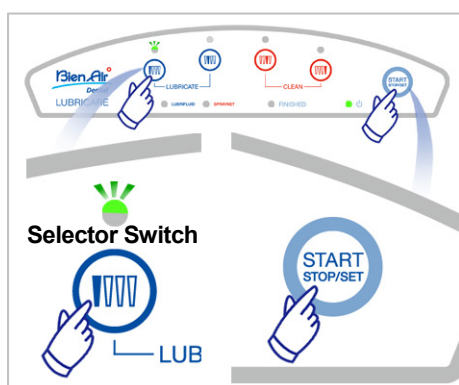


2. Hold down the Start Switch (Start / Stop / Set Switch) and then hold down a Selector Switch for 2 seconds. Release it when a beep sounds. The Selector Lamp will blink to show that you can now change its setting.

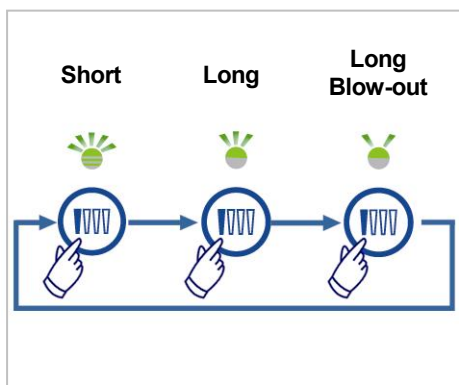

3. Now press the Selector Switch to change modes. Each press of the Selector Switch will change the modes in sequence.

* The lamp's blinking rate will change.

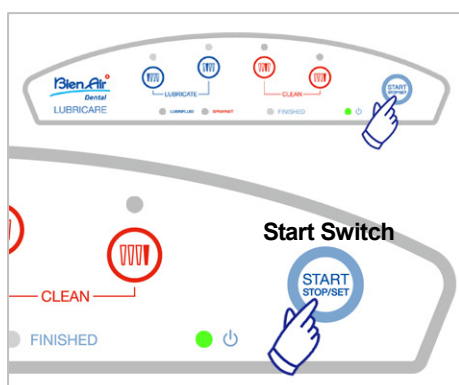

4. Check the time setting by how fast the lamp blinks, and then press the Start Switch for 2 seconds. A double beep signals that the change has been made. After this the lamp goes out.

Usage Note


♦ *Each Select Switch must be set separately using the above procedure.*

Coupling Replacement

1. Loosen the coupling's nut and take it off.



2. Fit the replacement coupling onto the nozzle and tighten up its nut.

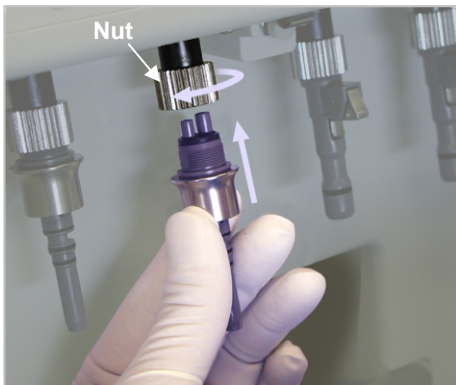

Usage Note

- ♦ *Make sure the nut is securely finger tight. Otherwise oil and air will escape and the maintenance will not be properly performed.*
- ♦ *Do not tighten the nut with pliers or other tools. This will damage the coupling.*

⚠ CAUTION

For handpieces not made by Bien-Air Dental SA, do not fail to use the correct LUBRICARE coupling or the coupling for the main tube.

Do not use couplings made for maintenance equipment other than LUBRICARE. This could result in release of oil mist or poor handpiece maintenance.

Replacing Fuses




The power supply code is detached from the main body.

The fuse holder is pulled out.

The fuse is exchanged for the new article.

Fuse: 250 V 2 A Slow-Blow and High breaking Type

The fuse holder is inserted.

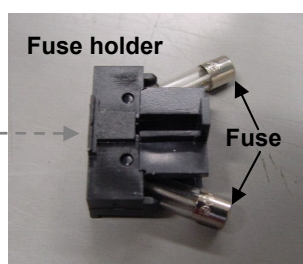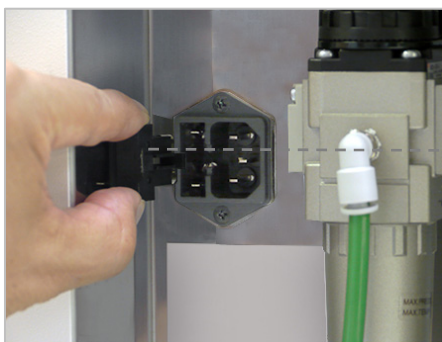
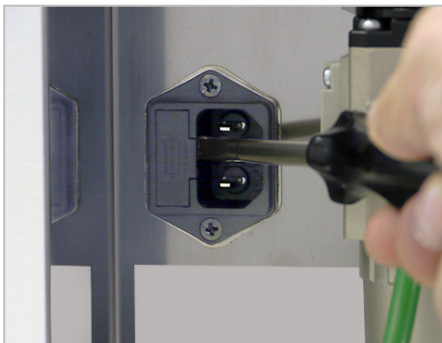


6. Maintenance and Replacement Parts

WARNING

Wear surgical gloves to operate and clean LUBRICARE.

Remove and Clean Front Door

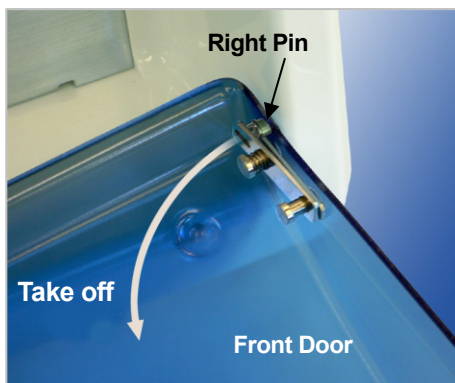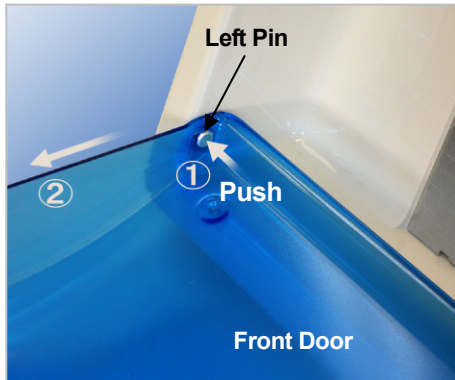

1. Push the pin on the left side in the direction of the arrow in the photo.
Push the pin all the way in and then pull the door forward.

2. Take the door off the pin on the right side.

CAUTION

Take care not to cut your finger on the hardware.

Usage Note

- ♦ *Hold the door securely and pull it off carefully. Excessive force could damage it.*

- * Fit the door back onto the pin on the right side to replace it.

3. After removing the door wash it with water and a neutral detergent or wipe it with Ethanol for Disinfection (Ethanol 70 vol% to 80 vol%).

Usage Note

- ♦ *Do not use hot water; this could damage the plastic.*
- ♦ *Use only soft gauze to wipe the door.*



4. Clean the main body with a neutral detergent or Ethanol for Disinfection (Ethanol 70 vol% to 80 vol%).

⚠ CAUTION

Take care not to cut your finger on the couplings or other parts.


Usage Note

- ◆ *Do not use disinfectants that contain chlorine-based cleaners or benzalkonium chloride or chlorinated aromatics. These solutions could discolor or damage plastic and corrode metal.*


Oil Pad Replacement

- * Inspect the oil pad once a week and replace it before oil reaches the top.

Remove the oil pad case.

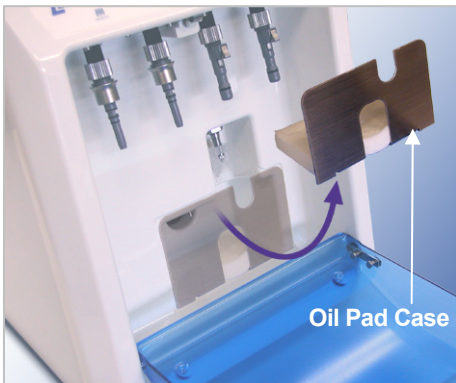

Oil Pad Case

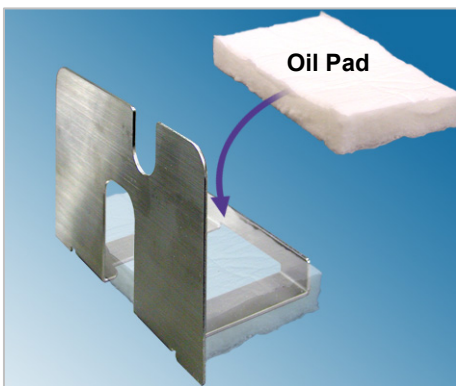

Oil Pad


Take the pad out of the case and replace it with a new one.

⚠ CAUTION

Dispose of used oil pads according to regulations for medical waste.

Usage Note


- ◆ *An oily film will form on LUBRICARE if the pad is not replaced.*
- ◆ *A single oil pad can absorb the oil from about 2 spray cans. However, the absorbed oil could start leaking from the pad if it is not replaced soon enough.*

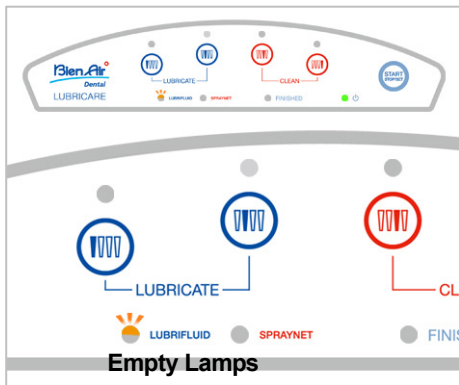

Spray Can Replacement

Replace a can if its empty lamp starts to blink.
(See pages 13 “(6) Spray Cans”.)

Empty Drainage from Air Regulator

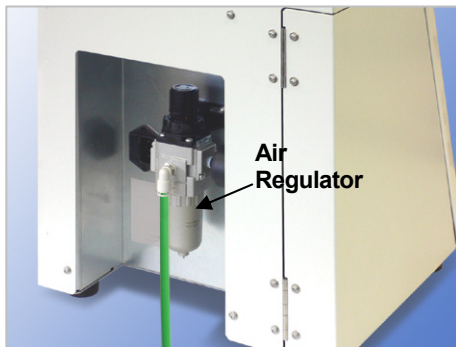


- * Empty the drainage from the air regulator about once a month.



1. Pull up the dial on the air regulator and turn it until the pressure drops to 0 MPa.

⚠ CAUTION

If the drain valve is opened without releasing the pressure inside the regulator, the water in the drain could spray out into your eyes.

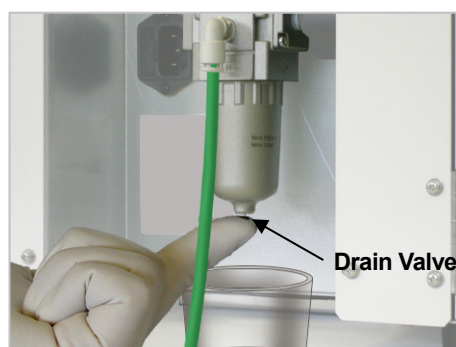

2. Put a container under the regulator and then press the drain valve.
The dirty water in the regulator will drain out.

⚠ CAUTION


If water is inadvertently mixed with the air, handpieces could be severely damaged.



3. After draining the regulator, turn the dial back to adjust the pressure and push it back down.

Cleaning the Tray


Pull out the tray and clean it with a neutral detergent or Ethanol for Disinfection (Ethanol 70 vol% to 80 vol%) about once a month.

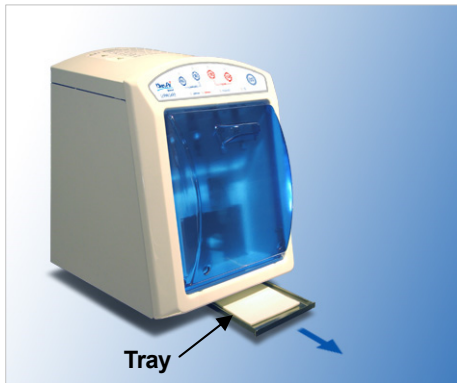

Usage Note

- ◆ *The tray cannot be autoclaved. Autoclaving will cause it to rust.*
- ◆ *Do not use disinfectants that contain chlorine-based cleaners or benzalkonium chloride or chlorinated aromatics. These solutions will corrode the tray.*
- ◆ *If this tray is not cleaned, an oily film will collect on the outside of LUBRICARE.*

Throw away the old filter and replace it with a new one.

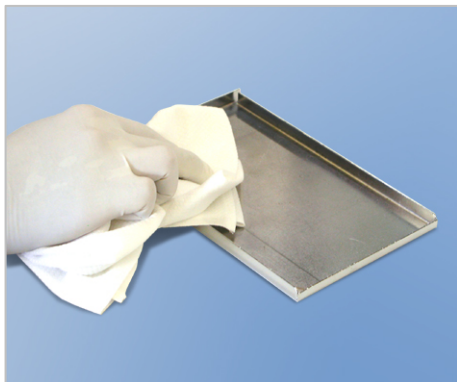

⚠ CAUTION


- *Dispose of used filters according to regulations for medical waste.*
- *Do not fail to push the tray all the way in. Otherwise, oil could leak out or oil mist might be released.*

Cleaning Couplings

Loosen its nut and take the coupling off. Clean it with Ethanol for Disinfection (Ethanol 70 vol% to 80 vol%).

Usage Note

- ♦ *Take care not to damage the O-ring.*
- ♦ *Couplings cannot be autoclaved. Autoclaving could crack the plastic parts.*

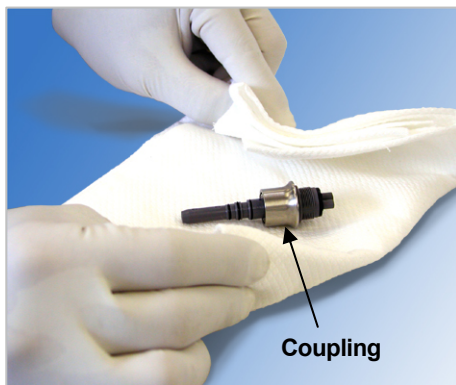

- * Refer to page 24 “Coupling Replacement” for how to attach couplings.

7. Maintenance and Inspection

Regular Inspection

- * Maintenance and inspection are generally considered to be the duty and obligation of the user, but if, for some reason, the user is unable to carry out these duties, he may rely on a qualified medical device serviceman. Contact your local dealer or Bien-Air Dental SA for details.

1. Power Supply Cord
Inspect visually for wear and broken wires.
2. Main and Operation Switches
Turn main switch on and check that main lamp lights up. Check that the unit operates correctly.
3. Oil Absorption Pad
Open case and replace oil pad if oil has seeped all the way to its top. (Replace after using 2 cans of spray.)
4. Drain Air Regulator
Check and drain air regulator. (once a month)
5. Tray Filter
Replace tray filter. (once a month)

Regularly Required Replacement Parts

Fuse

Other Replacement Parts

Oil Absorption Pad

Tray Filter

8. Troubleshooting

Use the check list below if the apparatus does not seem to be working properly.

- * Contact your local dealer or Bien-Air Dental SA, if the apparatus does not work normally even after inspection, adjustment or parts replacement, or if the dentist cannot perform maintenance procedures himself.
- * Before inspection and adjustment, make sure the power is on.

Problem	Cause	Solution
Main power lamp does not light up.	Power supply cord not properly plugged into supply socket.	Plug in properly.
	Power supply cord not properly plugged into LUBRICARE.	
Does not start when Start Switch is pressed. (1) Double beep sounds. (2) Spray can Empty Lamp is blinking.	Front door may be open.	Close front door.
	Empty can	Replace can.
	Can is not properly installed.	Fix can installation.
	Compressor air not connected.	Connect air.
	Air regulator not set right.	Set regulator properly.
Unusual odor. Surfaces are unusually dirty.	Use without connecting handpiece.	Select only couplings that have a handpiece connected.
	Oil pad is saturated.	Clean LUBRICARE and replace oil pad.
	Tray filter is saturated.	Clean LUBRICARE and replace filter.
	Spray can not set right.	Install can properly.
	Air regulator not set right.	Set regulator properly.
Cannot install spray can. Oil leaks from spray can.	Bent spray can nozzle pin.	Replace can.
	Wrong stand for spray can.	Use proper stand for can.
	Spray can by maker other than Bien-Air Dental SA.	Use a Bien-Air Dental SA spray can.
Oil mist is released. Equipment surfaces get oily.	Handpiece type	Depending on the type of handpiece and coupling, there could be a release of oil mist. In this case, lower the pressure to 0.3 MPa.
	Air pressure is too high.	Reduce air pressure, but not less than 0.3 MPa.
	No burr inserted in straight attachment.	Refer to the user instructions for the handpiece. * For some handpieces, an oil mist may be released if a burr is not inserted.

9. Service Contacts

* For repair or other types of service contact your local dealer or Bien-Air Dental SA.

* Disposal of Medical Devices

Any medical devices which could possibly be contaminated must be first decontaminated by the responsible doctor or medical institution and then be disposed by an agent licensed and qualified to handle medical and industrial waste.



J. MORITA MFG. CORP.

680 Higashihama Minami-cho, Fushimi-ku, Kyoto 612-8533, Japan

EU Authorized Representative under the European Directive 93/42/EEC



MEDICAL TECHNOLOGY PROMEDT CONSULTING GmbH

Altenhofstraße 80, 66386 St. Ingbert, Germany
T +49. 6894 581020, F +49. 6894 581021

The authority granted to the authorized representative, MEDICAL TECHNOLOGY PROMEDT Consulting GmbH, by J. MORITA MFG. CORP. is solely limited to the work of the authorized representative with the requirements of the European Directive 93/42/EEC for product registration and incident report.

10. Appendix – Electromagnetic Declaration

This device conforms to IEC 60601-1-2: 2007, the relevant international standard for electromagnetic compatibility (EMC).


The following is the “Guidance and Manufacturer’s Declaration” which is required by IEC 60601-1-2: 2007, the relevant international standard for electromagnetic compatibility.

WARNING

- ***LUBRICARE (hereafter referred to as the HIM-1-BA) needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the ACCOMPANYING DOCUMENTS.***
- ***Portable and mobile RF communications equipment can affect the HIM-1-BA.***
- ***Use of parts other than those accompanied or specified by manufacturer or distributor may result in increased EMC emissions or decreased EMC immunity of the HIM-1-BA.***
- ***The HIM-1-BA should not be used adjacent to with other equipment. If adjacent use is necessary, the HIM-1-BA should be observed to verify normal operation in the configuration in which it will be used.***

Guidance and Manufacturer’s Declaration – Electromagnetic Emissions		
The HIM-1-BA is intended for use in the electromagnetic environment specified below. The customer or the user of the HIM-1-BA should assure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment – Guidance
RF emissions CISPR 11	Group 1	The HIM-1-BA uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The HIM-1-BA is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
*Harmonic emissions IEC61000-3-2	Class A	
*Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

Guidance and Manufacturer's Declaration – Electromagnetic Immunity			
The HIM-1-BA is intended for use in the electromagnetic environment specified below. The customer or the user of the HIM-1-BA should assure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±2, 4, 6 kV contact ±2, 4, 8 kV air	±2, 4, 6 kV contact ±2, 4, 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
*Electrical fast transients/bursts IEC 61000-4-4	±2 kV for power supply lines	±2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
*Surge IEC 61000-4-5	±0.5, 1 kV line(s) to line(s) ±0.5, 1,2 kV line(s) to earth	±0.5, 1 kV line(s) to line(s) ±0.5, 1,2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
*Voltage dips, short interruptions and voltage variations on power supply lines IEC 61000-4-11	<5% U_T (>95% dip in U_T) for 0.5 cycle 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles <5% U_T (>95% dip in U_T) for 250 cycles	0% U_T (>95% dip in U_T) for 0.5 cycle 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles 0% U_T (>95% dip in U_T) for 250 cycles	Mains power quality should be that of a typical commercial or hospital environment. If user of the HIM-1-BA requires continued operation during power mains interruptions, it is recommended that the HIM-1-BA be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3.15 A/m	3.15 A/m	Power frequency magnetic field should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note U_T is the a.c. mains voltage prior to application of the test level.			

Guidance and Manufacturer’s Declaration – Electromagnetic Immunity			
The HIM-1-BA is intended for use in the electromagnetic environment specified below. The customer or the user of the HIM-1-BA should assure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment – Guidance
Conducted RF IEC 61000-4-6	3.15 Vrms 150 kHz – 80 MHz	3.15 Vrms	<p>Portable and mobile RF communications equipment should be used no closer to any part of the HIM-1-BA, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance</p> $d = \frac{3.5}{3.15} \sqrt{P}$ $d = \frac{3.5}{3.7} \sqrt{P} \quad 80 \text{ MHz} - 800 \text{ MHz}$ $d = \frac{7}{3.7} \sqrt{P} \quad 800\text{MHz} - 2.5 \text{ GHz}$ <p>Where <i>P</i> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <i>d</i> is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range. ^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
Radiated RF IEC 61000-4-3	3.7 V/m 80 MHz – 2.5 GHz	3.7 V/m	
NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.			
NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicated theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the HIM-1-BA is used exceeds the applicable RF compliance level above, the HIM-1-BA should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the HIM-1-BA .			
^b Over the frequency range 150 kHz to 80MHz, field strengths should be less than 3.7 V/m.			

Recommended Separation Distances between Portable and Mobile RF Communications Equipment and the HIM-1-BA.

The **HIM-1-BA** is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the **HIM-1-BA** can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the **HIM-1-BA** as recommended below, according to the maximum output power of the communications equipment.

Rated Maximum Output Power of Transmitter (W)	Separation Distance According to Frequency of Transmitter (m)		
	150 kHz – 80 MHz $d = 1.11 \sqrt{p}$	80 MHz – 800 MHz $d = 0.95 \sqrt{p}$	800 MHz – 2.5 GHz $d = 1.89 \sqrt{p}$
0.01	0.11	0.10	0.19
0.1	0.35	0.30	0.60
1	1.11	0.95	1.89
10	3.51	3.00	5.98
100	11.11	9.50	18.9

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

*Basic Performance (Criteria):

- Noise does not cause the HIM-1-BA to malfunction.

*Cable Length:

- AC Cable: 2.5 m (2-Wire), 3.0 m (3-Wire)

BIEN-AIR ADDRESSES

Bien-Air Dental SA
Länggasse 60
PO Box
2500 Bienne 6, Switzerland
Tel. +41 (0)32 344 64 64
Fax +41 (0)32 344 64 91
dental@bienair.com

Bien-Air
Deutschland GmbH
Jechtlinger Straße 11
79111 Freiburg, Deutschland
Tel. +49 (0)761 45 57 40
Fax +49 (0)761 47 47 28
ba-d@bienair.com

Bien-Air España SA
Entença, 169 Bajos
08029 Barcelona, Spain
Tel. (+34) 934 25 30 40
Fax (+34) 934 23 98 60
ba-e@bienair.com

Bien-Air USA, Inc.
Medical Technologies
5 Corporate Park
Suite 160
Irvine, CA 92606, USA
Tel: 1-800-433-BIEN
Tel: 949-477-6050
Fax 949-477-6051
ba-usa@bienair.com

Bien-Air France Sàrl
19-21, rue du 8 Mai 1945
94110 Arcueil, France
Tel. +33 (0) 1 49 08 02 60
Fax. +33 (0) 1 46 64 86 68
ba-f@bienair.com

Bien-Air Italia s.r.l.
Via Vaina 3
20122 Milano, Italia
Tel. +39 (02) 58 32 12 51/52/54
Fax +39 (02) 58 32 12 53
ba-i@bienair.com

Bien-Air UK Limited
Arundel House
Unit 1 - Ground Floor
Amberley Court, Whitworth Road
Crawley, West Sussex, RH11 7XL
England
Tel. +44 (0)1293 550 200
Fax+44 (0)1293 520 481
ba-uk@bienair.com

Beijing Bien-Air
Medical Instrument
Technology Service Co. Ltd.
Room 907, The Exchange Beijing,
No 118 Jian Guo Lu Yi,
Chao Yang District,
Beijing 100022, China

北京彼岸医疗器械
技术服务有限公司
北京市朝阳区建国路
乙 118 号招商局中心
京汇大厦 907 室
电话: +86 10 6567 0651
传真: +86 10 6567 8047
ba-beijing@bienair.com

www.bienair.com